

X-RAY GENERATORS

CMP 200®

(VZW2556 Series)



CMP 200[®] series

Radiographic X-Ray generators that deliver the best price / performance ratio in the market".

High frequency performance in excess of 400kHz.


3 different models to choose from:

- 32 kW and 40 kW
 - 1 Phase 208/230 VAC (125kV) only
 - 3 Phase 208/230 VAC (125 kV) only
 - 3 Phase 400/480 (125 kV) standard (150 kV) optional
- 50 kW
 - 3 phase 208/230 VAC (150 kV) standard
 - 3 Phase 400/480 VAC (150 kV) standard

Cost-effective for all radiographic applications (hospitals, medical clinics, chiropractic and veterinary applications)

Small. Compact. Lightweight.
A cost effective choice.

X-RAY
**Generators
of Choice**

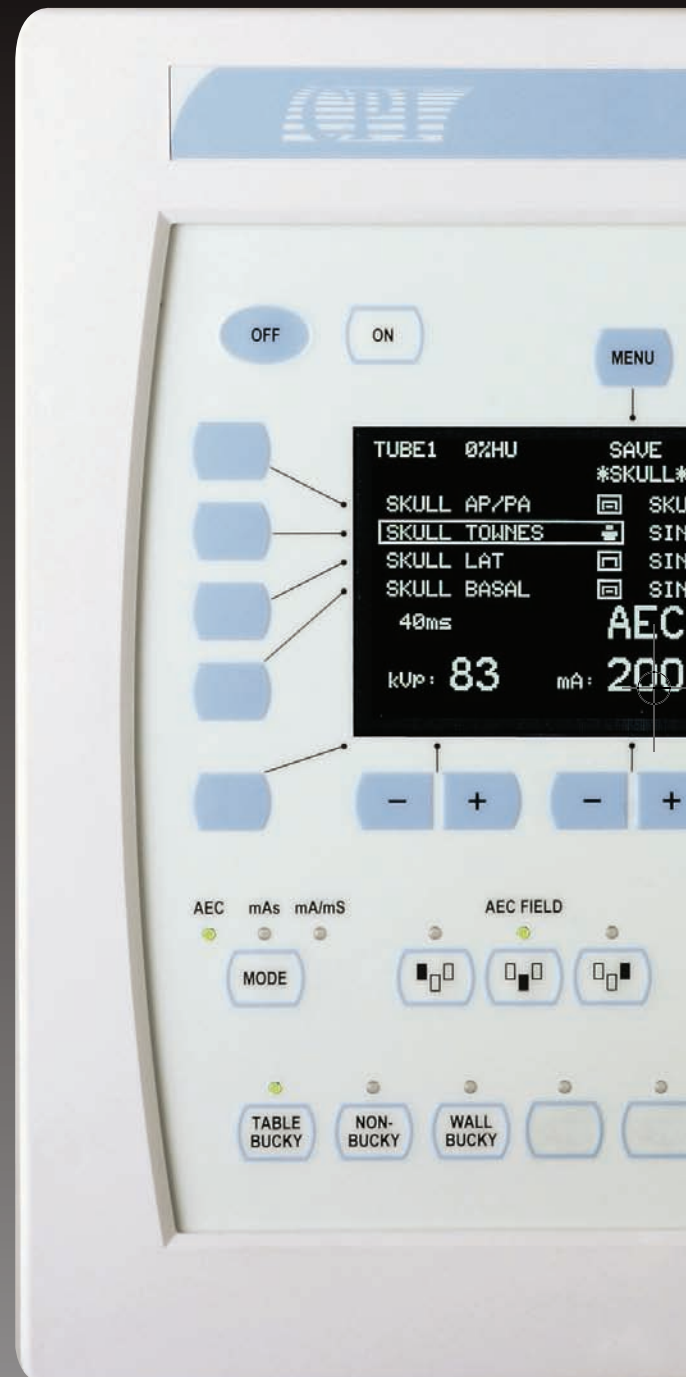



Look at these Features

- State-of-the-art 400kHz + FET-based inverter technology from the Indico 100® is also present in the CMP 200®.
- Anatomical programs available in several languages – and easily edited using the Operator Console
- Technique Selection:
 1. kV with AEC (Optional)
 2. kV/mAs (Standard)
 3. kV/mA/Time (Standard)
 4. kV/cm Thickness (Standard)
- Low kV ripple -- maximizes image quality by providing constant potential performance.
- CPI's output frequencies in excess of 400 kHz are the highest in the industry and produce the lowest kV ripple, resulting in the highest potential image quality and the lowest stored energy. This also means you get maximum X-ray tube protection and longest tube life.
- Operator Control Console available with international symbols or text.
- Wall Mount (Optional)
- Automatic Exposure Control (AEC) (Optional)
- Dose Area Product (DAP) Interface (Optional)
- Automatic Tube Calibration (Standard)
- Compact, cost effective design saves you money in transportation and installation.
- Hand Switch (Optional)
- Collimator Power Supply
- Tomography Interface (Optional)
- Multiple Comm. ports for easy interfacing



Hand Switch
(Optional)

World Class Product Support

Our greatest strength is our world class reputation for customer support. A dedicated support team is only a phone call away at 905-877-0161 – or at our North American toll-free number 1-888-CPI-XRAY (274-9729)

Call us. We are here to help you with all of your installation needs.



Communications & Power Industries Canada Inc.
communications & medical products division
45 River Drive / Georgetown, Ontario / Canada L7G 2J4
Hot Line Telephone: 1-888-CPI-XRAY
TEL: 905-877-0161 / FAX: 905-877-5327
E-MAIL: marketing@cmp.cpii.com
WEB: www.cpii.com



Some features shown may be optional. Please refer to Data Sheet.